

MiniScan® EZ

for...

FIBERGLASS
INSULATION



The color of fiberglass insulation represents the binding agent used by the insulation manufacturer. More importantly, the color is the manufacturer's personalized stamp on their product. So, insuring the color is exact to the manufacturers specification is brand critical. And for color measurement on-the-go, there is no better choice than the portability of the MiniScan® EZ 4500. The MiniScan® EZ 4500 has a large, easy-to-read display, is light weight with a rubberized handle for operator comfort, and ergonomically designed for one-handed operation and thumb-tip navigation of display functions. The onboard software includes most major industrial reflectance color scales and indices, stores 100 standards and 800 sample measurements, and can average up to 20 sample readings. MiniScan® EZ 4500 allows fiberglass insulation manufactures to easily and accurately evaluate the quality of their proucts. The MiniScan® EZ 4500 is compatible with HunterLab's EasyMatch® QC software and can be connected to a printer if hard copy data is desired.



**Color measurement
on-the-go**



**Easy one-handed operation
reduces operator fatigue**



**45°/0° geometry measures color the
way the eye sees color**



**Large, easy-to-read LCD
graphical display**



**Onboard software includes
most commonly used color
scales and indices**



measure
color
the way your
eye sees it

The Best Solution for Color Measurement of Fiberglass Insulation




Small and rugged – designed for portability

Providing operator ease-of-use, the lightweight MiniScan® EZ 4000 is ergonomically designed for one-handed operation, with thumb-tip navigation of all display functions, reducing operator fatigue and preventing sample reading errors. The MiniScan® EZ 4000 even provides a rubberized handle for increased grip and comfort. Its LCD display screen can be configured with four separate rotation choices: 0°, 90°, 180°, and 270°.



45°/0° measurement geometry

MiniScan® EZ 4500 measures in the visible spectrum (400nm-700nm) using 45°/0° geometry for the most accurate and consistent color measurement of fiberglass insulation.

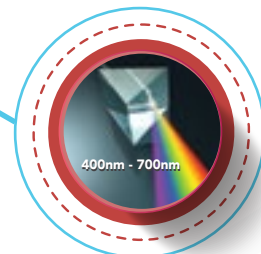


Onboard software includes everything needed for fiberglass insulation color measurements

The onboard MiniScan® EZ 4000 software includes all common color industrial scales and indices, CIE L*a*b*, CIE LCh, Hunter L,a,b and more. Store 100 product set-ups and up to 800 sample measurements including standards values and pass/fail tolerances. View color data, color difference data, spectral data, spectral data difference, spectral plot, and spectral difference plot on the large easy to read LCD graphical display.



Connect to a printer or download data PC software

Connecting to a USB printer allows hard copies of measurement data to be printed, or, download measurement data to HunterLab EasyMatch® QC PC software to collect, display and analyze your color data.

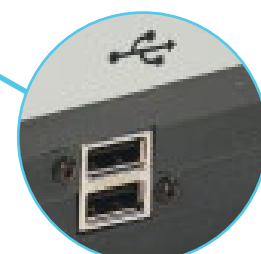

Keys for the successful measurement of Fiberglass Insulation



MiniScan® EZ 4500 comes ready to measure out of the box, and includes NiMH batteries, Battery charger, Calibrated instrument white tile, Certificate of traceability, Black glass, Green diagnostic tile, Dust cover, Carrying case, USB flash drive, and MiniScan® EZ users guide.

Optional Accessories

	D02-1014-427	45°/0°LAV Nose Cone with Screw on Polycarbonate Port
	A13-1014-294	USB Flexible Keyboard
	A13-1014-259	USB Thermal Printer
	EZMQC	EasyMatch® QC Software