




Agera®

Patent Pending

for...

KRAFT PAPER

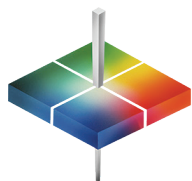

HunterLab

From matte to glossy paper finishes, the 45°/0° circumferential geometry is the best geometry for consistent and accurate color measurement results of both non-optically brightened and optically brightened Kraft Paper. Add to this the built in 60° gloss meter compliant with ASTM D523 and ISO2813, and you get a complete appearance result in a single measurement.

With modern data communications including USB, Ethernet, and Wi-Fi, and a capacitive high-resolution color touchscreen display, and easily to navigate embedded EasyMatch Essentials quality control software, Agera is the perfect solution for measuring the color of Kraft Paper.



0°/45°-c Circumferential Geometry and 60° Gloss in a Single Measurement



Full Spectrum LED Illumination with UV Control



On-Board Sample Viewing and Image Capture



Embedded EasyMatch Essentials Quality Control Software – No PC Required



Smart User Interface and Data Management – USB, Ethernet and Wi-Fi



measure
color
the way your
eye sees it

The Best Solution for Color Measurement of Kraft Paper




0°/45°c Circumferential Geometry and Built-in 60° Gloss Meter

45°/0° circumferential geometry is the best geometry for consistent and accurate color measurement results of Kraft paper because it replicates the human visual response to light and color. Add to this the built in 60° gloss meter compliant with ASTM D523 and ISO2813, and you get a complete appearance result in a single measurement

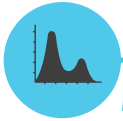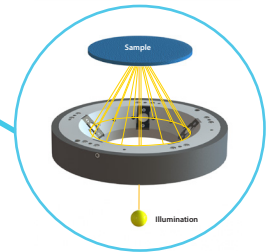


Full Spectrum, Balanced LED Illumination

Illuminating from 360-700nm, Agera allows you to measure and quantify fluorescent and optically brightened Kraft paper with simultaneous UV included and UV Excluded measurements.

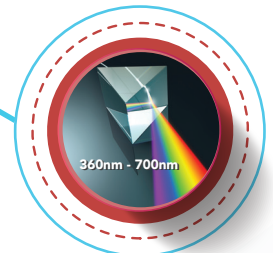

Sample Viewing and Image Capture

Agera Includes built-in 5-megapixel camera to view live on-screen images for precise sample placement and allows the image to be captured and stored as part of the data file.

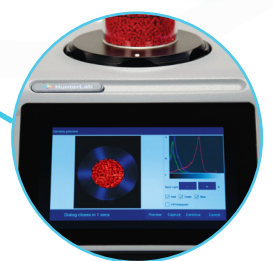

Embedded EasyMatch Essentials Quality Control Software


The onboard EasyMatch Essentials software is intuitive and easy to navigate. It offers applications to streamline production, allowing users to create their own workspaces to meet their specific application workflow needs. EasyMatch Essentials comes preloaded with most standard reflectance scales and indices widely used in industry.

Smart User Interface and Data Management - No PC Required

View measurement results and color graphics on the 7-inch 1280x800 capacitive industrial touchscreen, and easily communicate data where it needs to go via USB, USB on-the-go, Ethernet, Wi-Fi, and built in email.



Keys for the Successful Measurement of Kraft Paper

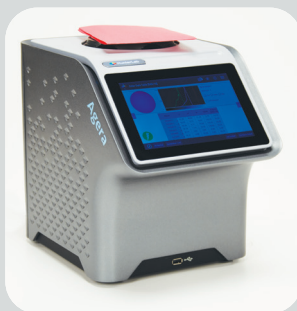

Agera includes everything needed to be successful measuring Kraft Paper in a port up configuration.

Optional Accessories



For Port-Forward sample measurements, the optional Sample Clamp should be ordered:



D02-1018-462 Sample Clamp